

Behind-The-Ear(BTE) Rechargeable Hearing Aid Manual

Please read the manual carefully to ensure proper use and keep the manual on file for future reference.


The Content of Component



- ① Rechargeable Case
- ② Hearing Aid
- ③ Ear Cap (large, medium, small)
- ④ Cleaning Brush
- ⑤ AA Alkaline Battery
- ⑥ AC Adapter


Maintenance Method

- Please clean your ear by using cotton swabs before you wear the hearing aid. If the hollow of sound is clogged by your earwax, the quality of volume will be inferior.
- Please clean out the earwax adhered on the hearing aid. The easiest location on which earwax to be stuffed is around the hollow of sound.

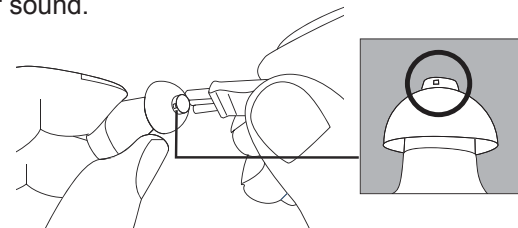


Diagram

Rechargeable Case

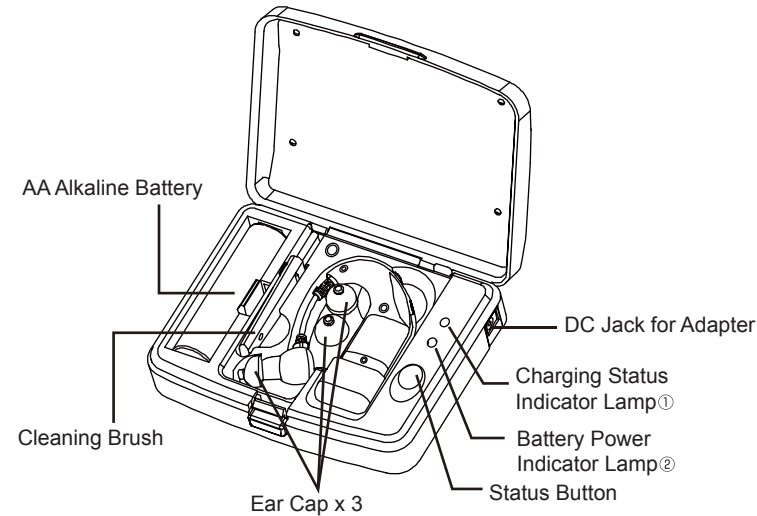


Hearing Aid

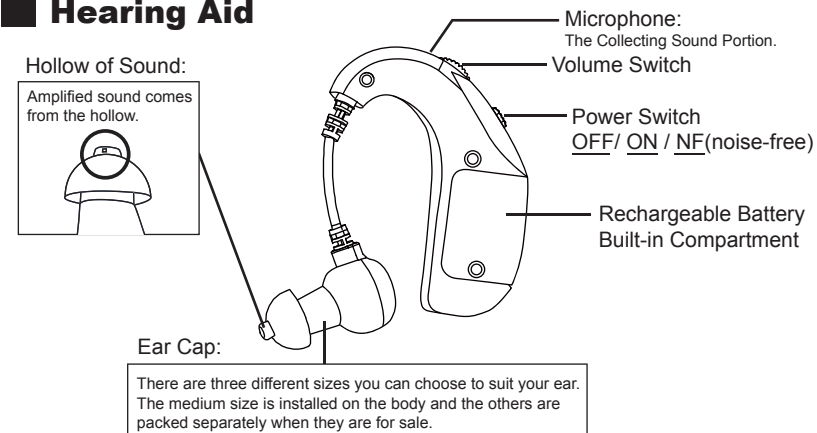

Setting Method of Hearing Aid

- **NF (Noise-Free) - Low frequency background noise excluded**
By state of the noisy environment, please set this function to hear in comfort.
- **ON** - When in nature and quiet environment
- **OFF** - Power off


Note: Before and After Each Use

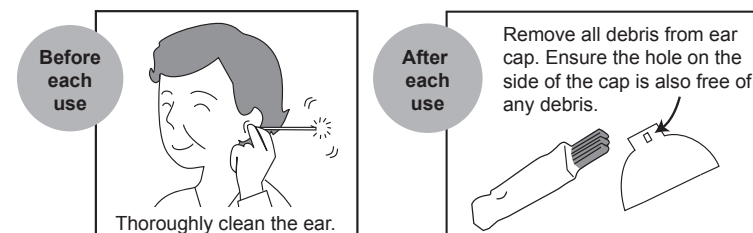

Recharge Method of Hearing Aid


■ The occasion of using Battery to recharge

- Install one attached battery in the compartment.
- Please observe whether the illuminant of the lamp ① is "GREEN" after you push on the status button.
- Please charge it for 7 hours at least.

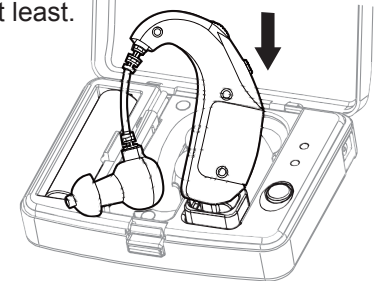


■ The occasion of using AC Adapter to recharge


- Insert the AC Adapter.
- Please observe whether the illuminant of the lamp ① is "GREEN" after you push on the status button.
- Please charge it for 5 hours at least.

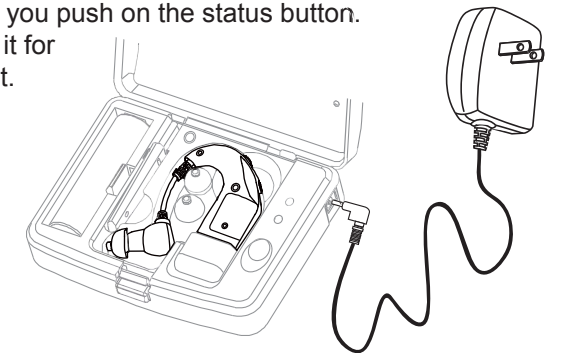


Rechargeable Lamp ① • ② (Illuminant will be shown after you press the status button)

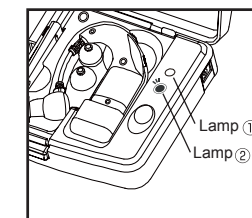


The occasion of "GREEN" illuminants come from Lamp ① and Lamp ② : It indicates under charge condition & stable power supply.

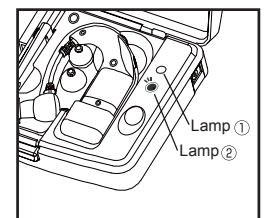

The occasion of "NO" illuminant comes from Lamp ① : The Hearing Aid is full charged.

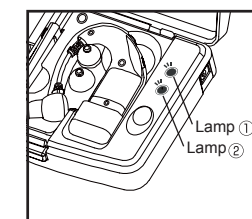

The occasion of "GREEN" illuminant comes from Lamp ② : It shows the stable power supply from the AA alkaline battery or AC adapter and the unit is full charged.

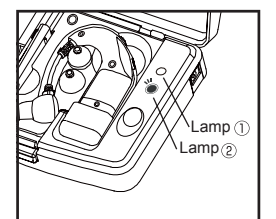

The occasion of "RED" illuminant comes from Lamp ② : It implies that either the AA alkaline battery is not properly installed or under low power. Please replace a new battery or use the AC adapter.

Before Request for a Repair Service

Please have an examination by yourself before sending it to our Repair Service Center. If all methods are ineffectual, please contact the center.

Can't hear any sound :

1. Does the power switch on?

Please set the power switch in the middle position to turn the power on.

2. Did the ear cap clog with earwax?

Please use our accessory - "cleaning brush" to clean up your ear cap.

3. Has the attached alkaline battery worn out completely?

Please replace with new alkaline battery and have the product charged.

4. Suggest the user to visit his/her physician to conduct a medical check periodically, if there is a reception problem.

5. Please fully charge this unit before the first use.



NOTICE

- * Please don't try other methods written out of the manual we enclosed. It may cause accidents or breakdowns.
- * Please don't let infants or people lacking the ability of self-expression. It may cause accidents or troubles.
- * Please don't use it while you are sleeping. It may cause accidents or physical discomforts.
- * Please don't wear it and do the fierce exercise at the same time. It may cause a breakdown or damage when you dropped it accidentally.
- * Please don't use this product together with your mobile phone or PHS. The sound will be lower or other disturbance caused because of the effect of noise resulted from radio wave.
- * Please don't moisten it .If you do, please have it inspected. Further, please don't wear it while taking shower or washing your hair. Additionally, please be cautious to keep it dry while washing your face or using hair spray. Or it may cause breakdowns.

Note: Concerning the disposal of hearing aids or batteries, please obey the regulations of the local administration.



Very Important Notice

~ Recharge Any Time~

For properly use under normal condition, it is able to be recharging & discharging for 500 times.

In order to extend the usage life of the rechargeable battery, here we strongly recommend users to charge the hearing aid whenever you don't use it.

Please note that there is no need to recharge your hearing aid until it is almost out of power. To charge your hearing aid whenever you don't use it will be helpful to maintain the rechargeable battery & extend the usage life for recharging & discharging **up to 500 times.**

According to the character of rechargeable battery, when it is almost run out of power, the sound from hearing aid will reduce. Please don't treat it as any defective, just put your hearing aid back to the charging case or use the AC adapter to recharge your hearing aid.

Specification & Function

General Title of Commodity : BTE Rechargeable Hearing Aid

Purpose & Function : For assisting people of hearing obstacle

Standard Frequency : 1200 ~ 1600Hz

Max. Output Sound Level Gain : 25 ± 5dB

Highest Sound Pressure Level for Input SPL 90dB : 120dB ± 10dB below 130dB

Equivalent Input Noise Level : Below 29dB

Built-in Rechargeable Battery Current : Below 5mA / 2.5V

Rechargeable Time : Up to 500 times in total

Continuous Use : 12-16 hours when full charged

Operating Temperature / Humidity : +10°C ~ +40°C, 30% ~ 85% RH

Storage Temperature / Humidity : -10°C ~ +60°C, 30% ~ 85% RH

Size : <Hearing Aid> 42(L)x14(W)x19(H) mm

<Rechargeable Case> 90(L)x70(W)x24(H) mm

Weight : <Hearing Aid> 6.8g <Rechargeable Case> 67g

Material : <Rechargeable Case> ABS <Ear Cap> silicone <AC Adapter> ABS

In The Situation Like


1. Having Conversations:

Keeping the average volume is fine.




2. Watching TV :

Let's enjoy your favorite program! You won't interfere with family at all.

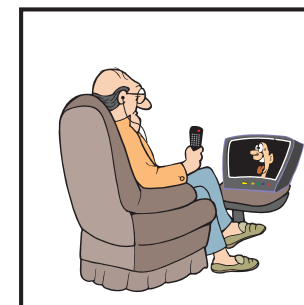

3. Bird-watch:

Clearly hear "Chirrup" in the distance.

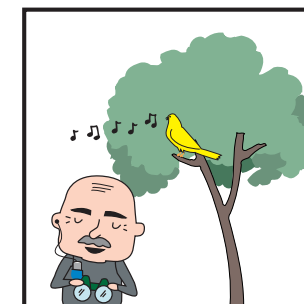

4. Listening to Speech:

Even sitting behind, you have no trouble listening a whole speaking.

